




# Terrafors

A solid state, in situ sterilisable bioreactor

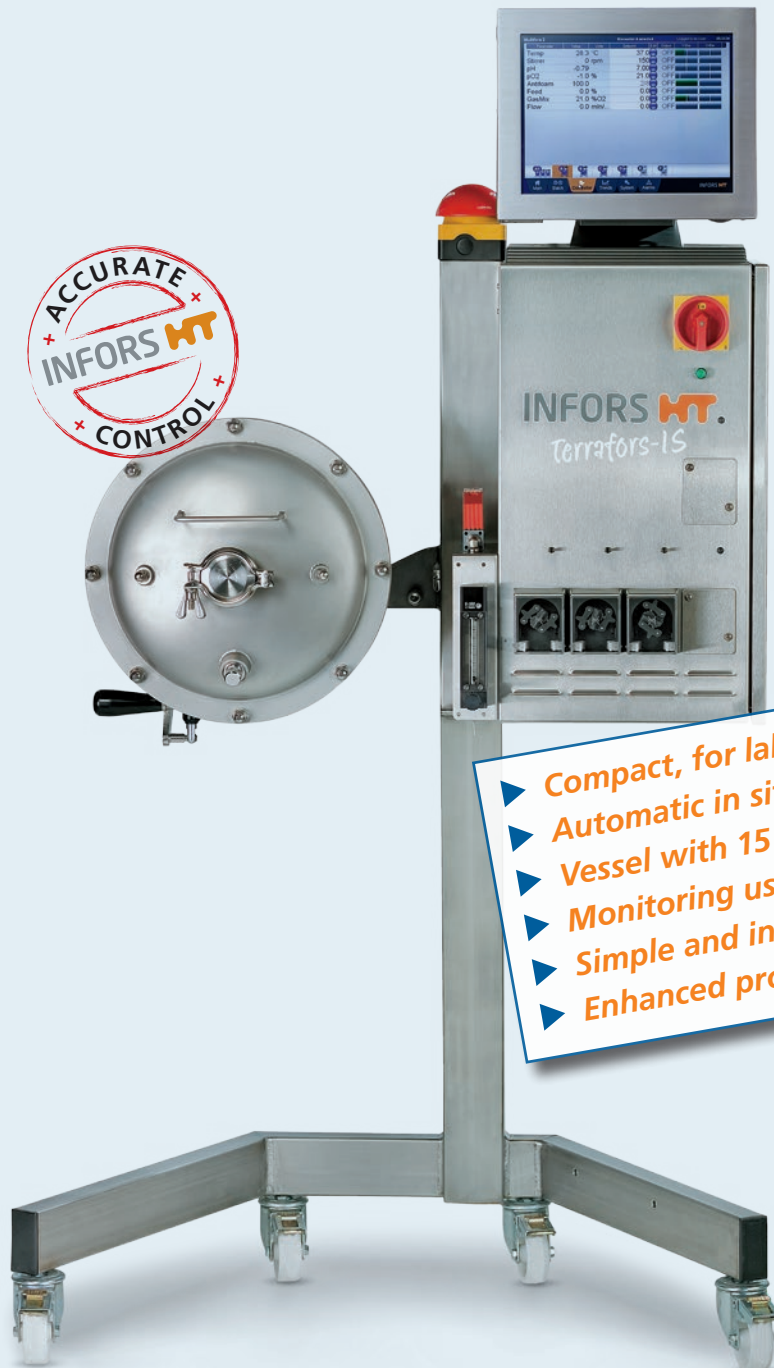

- ▶ Compact, for laboratory use
- ▶ Automatic in situ sterilisation
- ▶ Vessel with 15 L total volume
- ▶ Monitoring using exit gas analysis
- ▶ Simple and intuitive local control
- ▶ Enhanced process understanding



[www.infors-ht.com](http://www.infors-ht.com)

May be subject to technical amendments

INFORS HT

# Powerful Solid State Bioreactor

## ▶ Compact, for laboratory use

Terrafors is mounted on a single column with the vessel to one side and the instrumentation to the other for minimum use of space in the laboratory.

## ▶ Automatic in situ sterilisation system

The automatic in situ sterilisation system uses direct steam injection with a maximum temperature of 130°C. An optional steam generator can be mounted on the base for ease of use.

## ▶ Vessel with 15 L total volume

The Terrafors vessel is made from 316L stainless steel and provides a total capacity of 15 L. Maximal working volumes include 3–4 kg of solids (e.g., soil), semi-solids (e.g., sludge) or up to 7 L of liquids (e.g., oil). Two ports on the top plate are included for substrate addition/removal, and sampling.

## ▶ Monitoring using exit gas analysis

With the additional INFORS HT Gas Analyser the evolution of CO<sub>2</sub> can be related to rates of substrate utilisation and parameters such as respiratory quotient determined using oxygen and carbon dioxide ratios. These can, in turn, be used to e.g. control feed rates.

## ▶ Simple and intuitive local control

Terrafors is equipped with a high-performance touch screen controller. The responsive touch screen features a simple and intuitive local operation and allows easy access to all functions. Peripherals such as the INFORS HT Gas Analyser for exit gas measurement and analysis can be readily connected to the system.

**Up to 24 parameters:** Rotation speed/direction; Temperature (including in situ sterilisation); Feed; Gas flow; + several free channels, e.g., for exit gas analysis (O<sub>2</sub> and CO<sub>2</sub>) using the INFORS HT Gas Analyser or for pressure measurement and control up to 2 bar.

## ▶ Enhanced process understanding

Terrafors can be linked to the eve® software to plan, control and analyse your bioprocesses. The eve® software integrates workflows, devices, bio-process information and big data in a web-based platform that lets you organise your projects – no matter how complex they are.



### Applications

- **Production of enzymes and other biologically active compounds**
- **Bioremediation** of soils, slurries etc.
- **Biofuel production** through conversion of agricultural waste and similar material
- **Wastewater remediation**

### Key technical data

- Dimensions (W x D x H):** 990 x 970 x 1650 mm
- Total volume:** 15 L
- Working volumes:** 3–4 kg solids/semi-solids or 7 L liquids
- Speed range:** 0.1–10 rpm
- Temperature range:** ~5°C above coolant to 60°C (130°C for sterilisation)
- Standard parameters:** Rotation speed/direction, temperature (including in situ sterilisation), feed, gas flow
- Optional parameters:** Exit gas analysis, pressure control

Infors AG  
Headoffice, Switzerland

Rittergasse 27  
CH-4103 Bottmingen  
T +41 (0)61 425 77 00  
F +41 (0)61 425 77 01  
info@infors-ht.com

For more information and your local sales office please visit:

[www.infors-ht.com](http://www.infors-ht.com)

INFORS HT